


FlexX is our modern, ultra-high speed approach to ligand docking. Screen whole compound libraries within minutes to identify promising binders.

How does FlexX work?

The FlexX docking functionality in SeeSAR places a ligand into a binding site. It is based on an incremental construction algorithm. Ligands are split into so-called fragments, and an initial fragment (or combinations thereof) are placed into multiple places in the pocket – and scored using a simple yet very fast pre-scoring scheme. From the *n* solutions placed, the ligand is further built up, fragment by fragment, and the interim solutions are scored against each other. The best scored survive the process, and those are delivered to the user.

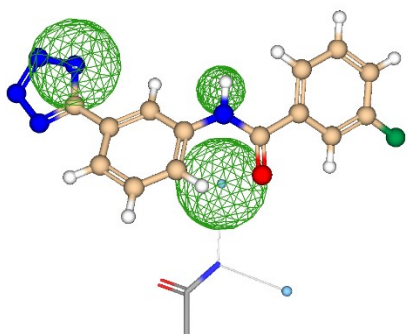


Advantages

- Docking for non-experts: no more receptor preparation – simply dock!
- Accurate binding mode prediction
- Dock gigantic libraries by using ultra-high speed docking (< 1s/ligand)

Virtual Screening

FlexX4 sets new records in virtual HTS. It runs in parallel and automatically uses all cores of your computer – thread-safe. You can screen a library of ~1,000 compounds in less than an hour on a laptop with 8 cores. FlexX-based docking led in many application studies to superior enrichment ratios [Warren 2006]. So you effectively find active compounds in large libraries. It uses ligand-based or protein-based pharmacophore constraints not only as a post-filter but as guidance to find in a vast solution space the ones which obey the constraints.

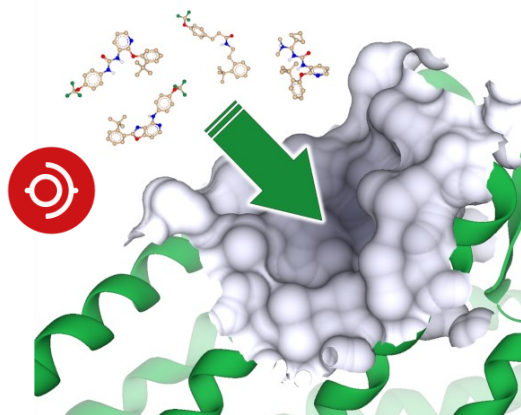

Accurate binding mode prediction.

For a protein with known three-dimensional structure and a small ligand molecule, FlexX4 accurately predicts the geometry of the protein-ligand complex within a few seconds. It has been validated, benchmarked and successfully used thousands of times as more than **10.000 citations** of the original scientific publication [Rarey 1996] impressively demonstrate.

Testimonials

- *“FlexX docking method is strait forward to automate. Everything is handeled by the tool itself – no preparation needed.”* – Hans Briem, Bayer, Germany.

Template-based docking.

Congeneric compound series can easily be docked. FlexX can superimpose a part of a ligand to dock onto a known crystal structure of a ligand. Template-based docking is performed by determining the Maximum Common Substructure (MCS) between the template and the compound-to-dock. Based on the MCS multiple overlays are generated to preserve the binding mode of the common core. The template need to have at least 5 heavy atoms.

Literature

A fast flexible docking method using incremental construction algorithm.

M. Rarey, B. Kramer, T. Lengauer, G. Klebe
1996 261 (3), 470-489
<https://doi.org/10.1006/jmbi.1996.0477>

A critical assessment of docking programs and scoring functions

G.L. Warren, C.W. Andrews, A.M. Capelli, B. Clarke, J. LaLonde, M.H. Lambert, M. Lindvall, N. Nevins, S.F. Semus, S. Senger, G. Tedesco, I.D. Wall, J.M. Woolven, C.E. Peishoff, M.S. Head
2006 49 (20), 5912-5931
<https://doi.org/10.1021/jm050362n>

Ultrafast de novo docking combining pharmacophores and combinatorics

M. Gastreich, M. Lilienthal, H. Briem, H. Claussen
2006 20 (12), 717-734
<https://doi.org/10.1007/s10822-006-9091-x>

Success Stories of Computer-Aided Design
 H. Kubinyi

2006 In book: Computer Applications in Pharmaceutical Research and Development, 377-424
<http://dx.doi.org/10.1002/0470037237.ch16>

